



Torso cushion

A comfortable support, e.g., during undressing for visits to the toilet.

Raising handles

Handles at two different heights that can also be adjusted for perfect adaptation to the height of each individual patient.

Seat frame

Provides stability in all direction.

Padded rear supports

Easy to fold up and shift aside, e.g., to facilitate undressing.

Safety lock

The rear support frame is easily secured in place with the reliable safety lock.

Lower-leg support

Height-adjustable and adapts automatically to the angle of the patient's legs.

Measurement scales

For easy adjustment of the lower-leg support, rear support frame and raising handles.

Width-adjustable base

Easily adjusted for wheelchair and toilet access and other situations.

Lockable wheels

The RollOn's wheels roll very easily. The rear wheels lock and unlock at a touch.

Foot tray

Durable non-skid design, easy to keep clean.


Right where you belong

RollOn is an easy-to-use manual support for raising that is specially designed for safe transfers to and from a sitting position. The RollOn is intended for patients who have considerable difficulty walking, but have good arm strength and the necessary leg strength to support their full weight.

Grasping the raising handles, the patient pulls himself up to a standing position and is surrounded by a stable, secure support and has comfortable, firm padded support for the knees and seat.

Lifts with care



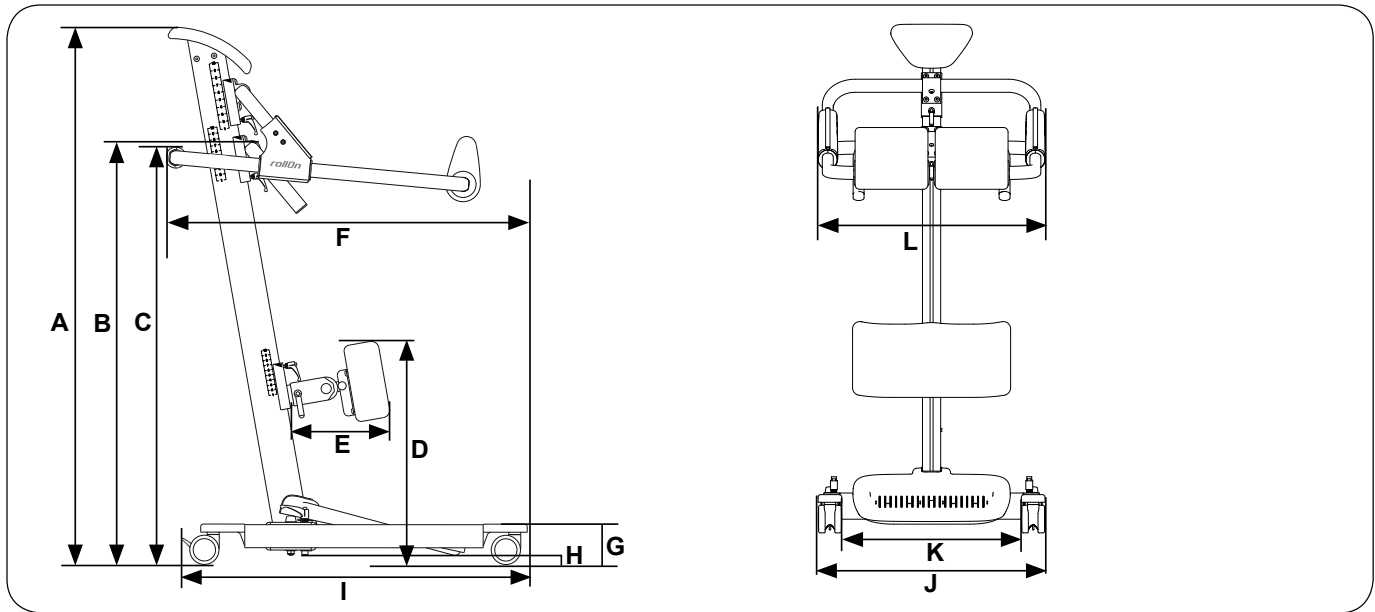

RollOn	Max. load	A		B*		C*		D*		E		F	G	H**	I	J		K		L	Weight total	Turning diameter
		Max	Min	Max	Min	Max	Min	Max	Min	Max	Min					Max	Min					
A	350	56	44	39	44	37	22	18	8	6	36	4	1	35	35	24	31	19	25	53	41	
B	160	141	112	100	112	95	57	46	20	16	92	11	2,6	89	90	62	78	49	63	24	105	

A: Max. load and weight in lbs. Measurements in inch. B: Max. load and weight in kg. Measurements in cm.
 * Max and min values within the measurement scale area. ** Floor clearance without load.

The RollOn is an excellent aid for transfers to and from the toilet or between bed and chair. It also makes it easier for the caregiver to assist a patient who has slid down in a wheelchair to achieve a better, more upright sitting position without the need for manual lifting. The RollOn contributes to active training in raising, without straining the caregiver.

To suit both short and tall patients, the RollOn is adjustable in a range of ways and there are several alternative hand-grips to facilitate raising. The padded rear supports can be easily folded up and aside, making it easy to get in and out of the RollOn and simplifying undressing and dressing during visits to the toilet.

Like many other Liko products, the RollOn is made of aluminium and is therefore very light and easy to keep clean.
 Max. load: 350 Ibs (160 kg)



RollOn complies with the Medical Device Directive for Class-1 products (MDD 93/42/EEC).

Made in Sweden

Liko is certified according to the ISO 9001 and ISO 14001 quality and environmental standards systems.



Lifts with care



www.liko.com

Manufacturer:
 Liko AB
 SE-975 92 Luleå
 Sweden
 info@liko.se

National Distributor in USA and Canada:
 Liko Inc • 842 Upper Union Street, Suite 4
 Franklin, MA 02038, USA
 Tel: 508-553-3993, 1-888-545 6671
 Fax: 508-528-6642
 info@likoinc.com

Liko Distributor in UK:
 Liko (UK) Ltd • Brunel Centre
 Stroudwater Business Park, Stonehouse,
 Gloucester, GL10 3RU
 Tel: +44 (0)1453 827272
 Fax: +44 (0)1453 828844
 info@liko.co.uk